

HotSwap™

by CytoSorbents

The new CytoSorb exchange system.

Effective hemoadsorption therapy often requires multiple adsorber exchanges, placing significant demands on patient care, treatment continuity and clinical staff.

In close collaboration with healthcare professionals, CytoSorbents has developed a pre-assembled blood line system that enables **seamless multiple adsorber exchanges.**



Features

- Integration into CRRT and Hemoperfusion (HP) platforms
- Designed and approved for CytoSorb exchange
- Bypass and adsorber lines for alternating blood flow
- Approved for 96h continuous use

i HotSwap can be suitable for all treatment settings when usage of more than one adsorber is indicated.

HotSwap can be suitable for clinical settings where

- CytoSorb dosing requires frequent changes
- CytoSorb timing requires early start
- Patient selection requires multiple usage



This might be applicable in the treatment of e.g.*



Hyperinflammation



Refractory shock




Rhabdomyolysis

Hyperbilirubinemia

*The final decision must be done at the discretion of the treating physician.

HotSwap is intended for a duration of use of up to 96 hours

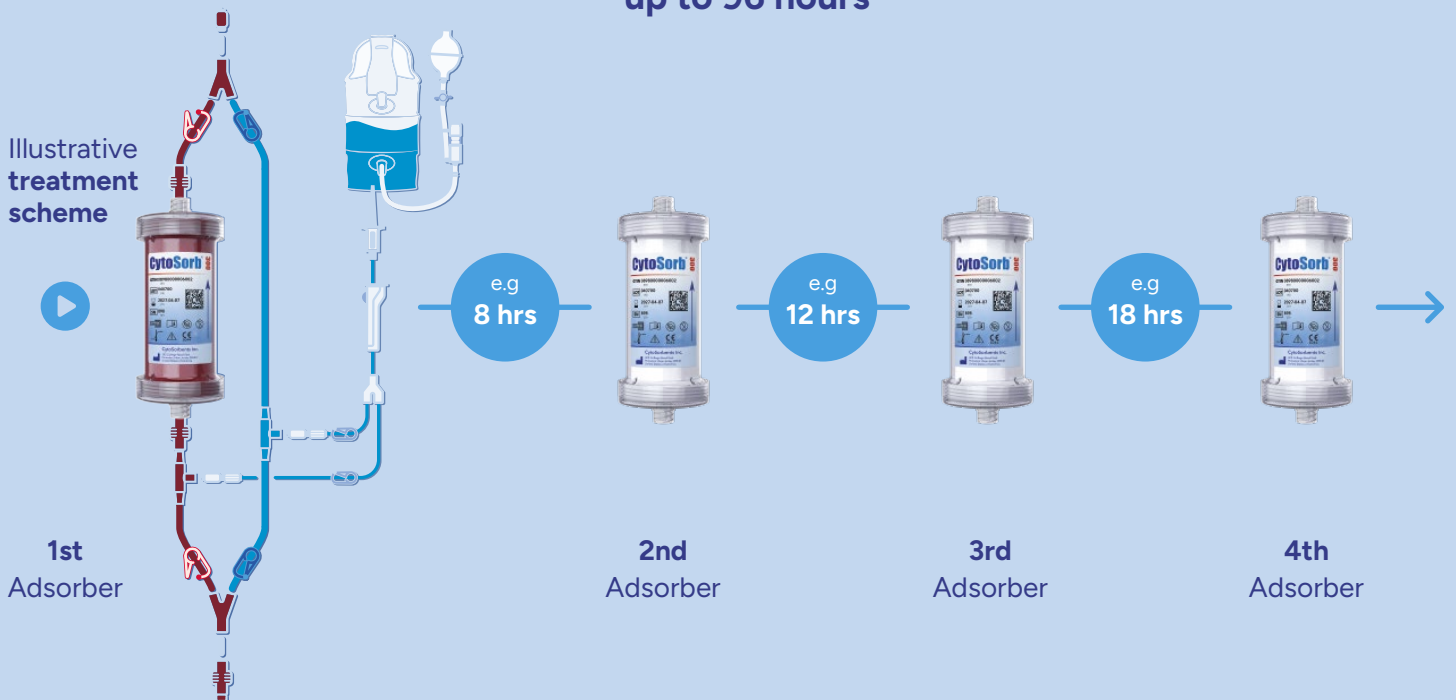

Key advantages of HotSwap for your daily practice



Safety

- No catheter manipulation
- Improved patient interface
- Simpler blood conservation



Flexibility

- Multiple adsorber changes
- Versatile Integration into extracorporeal circuits (EC)
- Flexible replacement scheme



Convenience

- Time saving workflow
- Intuitive handling steps
- Seamless circuit integration


Your advantages with HotSwap

- ✓ Convenient adsorber and bypass line combination
- ✓ Flexible adsorber exchanges with easy to handle steps
- ✓ Safe blood return with no circuit interruption


i The HotSwap system can be seamlessly integrated into the standard CRRT and extracorporeal circuits

For integration into CRRT, the lines shown from Adapter Sets 1a and 2a are required.


CRRT – Pre Filter

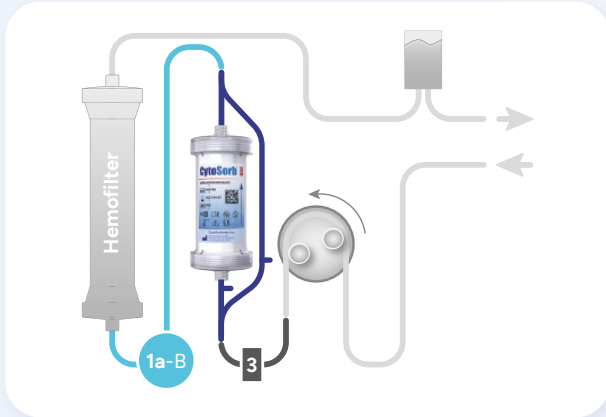

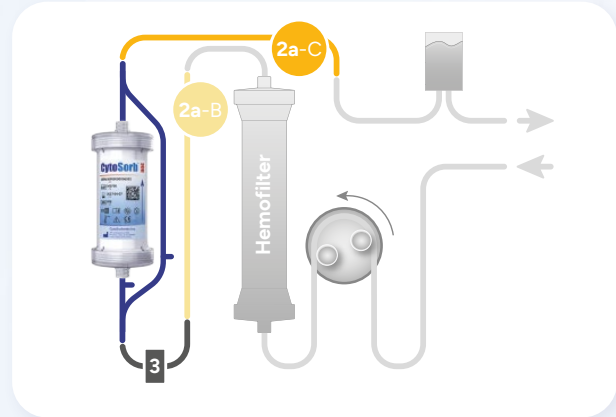


CRRT – Post Filter

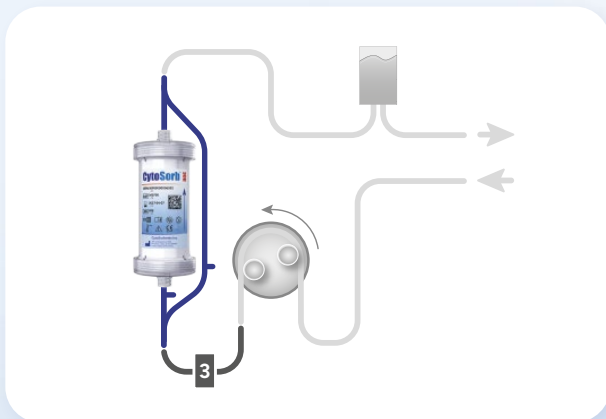

PuriFi / HP



- HotSwap
- Adapter Set 1a – B 
- Adapter Set 2a – B 
- Adapter Set 2a – C 
- 3 Bridge adapter

Visit us

For more information on how to apply HotSwap™, use our Quick Setup Guide:

→ www.cyto.zone/setup-hotswap/eng

